H3201

Maia Micro

Microdermabrasion / Microdermoabrasion / Microdermabrasion /
Microdermoabrasione / Mikrodermabrasion

EN	2
ES	5
FR	8
IT	11
DE	14

DDUUEETT EUROPE would like to thank you for choosing Sonia products. We hope that this instruction manual will help you get the most out of this product.

If you have any doubts regarding its use and/or operation, please contact your distributor.

BEFORE USE

- Please read the entire instruction manual carefully before using the product. Take note of any warnings included in this manual.
- Always keep this instruction manual handy and if you give the device to another person, hand over the manual as it is an essential part of the product.
- Some components of the product may differ from the image that appears in catalogues, websites, brochures, etc.
- In order to enhance the product, the manufacturer reserves the right to modify its features and/or components without prior notice.

CORRECT USE OF THE DEVICE

- This product has been designed to be used by beauty professionals at beauty salons and should not be used for other purposes.
- Please note that any misuse of the device entails the loss of the guarantee.
- Do not make any modification to the device without the authorisation of DDUUEETT EUROPE SL.
- DDUUEETT EUROPE SL shall not be held liable for any material or personal damage due to the misuse of the product.

CONTRAINDICATIONS

- This product should not be used on:
 - People with sensitive skin, lesions, inflammations, or any skin problem or disease.

PRECAUTIONS

- Before using it for the first time and after each use, check the device and the mains cable for any damage.
- Unplug the device from the mains before starting to clean it or undertaking any type of maintenance.
- Do not use the product outdoors or in places where it may be affected by humidity or inclement weather.
- If the power supply cable is defective, unplug the appliance immediately and inform your dealer or specialised technician.
- Position the power supply cable properly so that people cannot trip over it and make sure you plug it

into an earthed socket.

- At the end of the working day, disconnect the device and unplug the power supply cable from the mains.
- If you detect any anomaly in the operation of the product, please contact your distributor or specialised technician.

SAFETY INDICATIONS

- The device and accessories should be kept out of the reach of children.
- Only plug the device into a properly installed socket. The local mains voltage must correspond to the technical specifications of the device.
- The power socket must be accessible in case you have to quickly unplug the device at any given time.
- Always pull out the actual plug and not the mains cable.
- Never tamper with any electrical component or open up/repair any component of the device.
- Only your distributor or specialised technician should repair the device.
- Do not switch on the device if there is visible damage to it or to the mains cable.
- Never touch the device or the mains cable with wet hands.

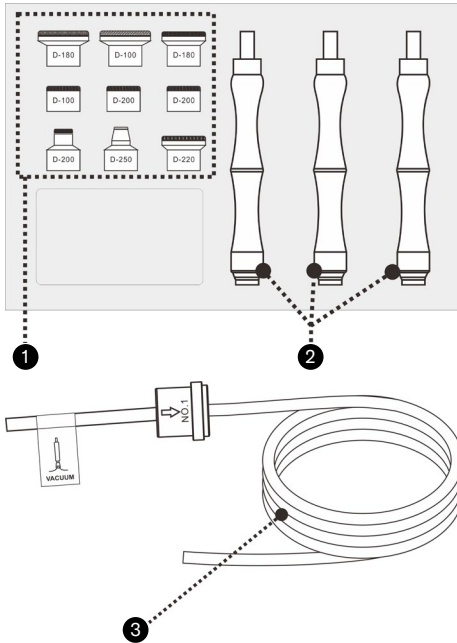
CONTENTS

- Main body x 1 unit
- Power cable x 1 unit
- Handpieces x 3 units
- Diamond-tipped heads x 9 units
- Flexible suction hose x 1 unit
- Sponge filters
- Instruction manual x 1 unit

TECHNICAL SPECIFICATIONS


Product	Microdermabrasion
Reference	H3201
Voltage	85-265 V
Input frequency	50/60 Hz
Power	50 W
Dimensions	28 x 17 x 11 cm
Net weight	1,80 Kg

COMPONENTS

- 1 Diamond-tipped heads
- 2 Handpieces
- 3 Flexible suction hose

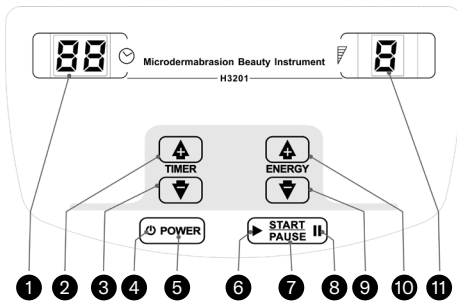
CONTROL PANEL

- 1 Timer
- 2 Increase time button
- 3 Decrease time button
- 4 Power on indicator
- 5 On / off button
- 6 Power on indicator
- 7 Start / Stop button
- 8 Pause indicator
- 9 Decrease power button
- 10 Increase power button
- 11 Power display

- 8 Pause indicator
- 9 Decrease power button
- 10 Increase power button
- 11 Power display

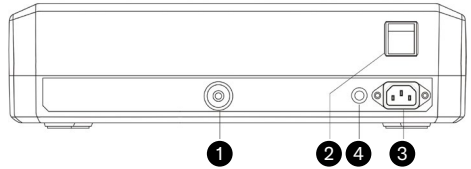
REAR PANEL

- 1 Connecting the suction hose
- 2 Power supply switch
- 3 Electrical power supply connector
- 4 Fuse

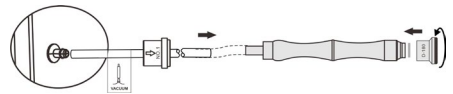CONNECTION AND ASSEMBLY INSTRUCTIONS

- First of all, remove two screws from the bottom of the main body. A sticker on the device indicates the need to remove them and highlights the two screws in question.
- Connect the power cable to the rear panel of the main body.
- Choose one of the diamond-tipped heads and connect it to the handpiece using the thread. Put the cotton filter inside the head first.
- Connect the suction hose to the corresponding outlet on the rear panel and on the end of the handpiece, as indicated in the picture. Note the direction of the arrow on the intermediate filter. It should be pointing towards the handpiece.

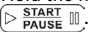
POWER AND CONFIGURATION

- Turn on the power switch on the rear of the machine and press **POWER** (power icon) on the control panel. The power indicator (power icon) will light up. The machine will be in standby mode and ready for operation. On the control panel:
 - The timer will show 15 minutes.
 - The power indicator will show 5.
 - On the **START PAUSE** (play/pause icon) button, the pause indicator (vertical bars icon) will light up.
- Set the timer using the **▲** or **▼** buttons to

choose the operating time of the machine, in a range of between 1 and 30 minutes.

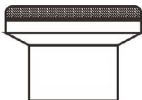
- Set the power using the  or  buttons to choose the power of the machine, in a range of between 1 and 9.
- Use the  button to start or stop the device.
 - When running, the  indicator will light up.
 - When paused or on standby, the  indicator will light up.

OPERATION


- Connect the device and leave it in standby mode with the  indicator lit.
- Set the power to the minimum and select the operating time.
- Hold the handpiece with your hand and press . The operating indicator  will light up and the device will start up.
- Place the device on the skin of the area to be treated to start the treatment, adjusting the suction power according to your needs.

Bear in mind:

- Clean the skin before starting the treatment.
- Make sure you have placed the cotton filter inside the selected head beforehand.
- To switch the heads or one handpiece for another, put the device in standby mode.
- The vacuum power should be kept under control, as a high power could cause skin lesions.
- Test first on your own hand, and then start with the patient once the vacuum power has been set.
- Regularly replace the cotton filter with a new one.
- The higher the number on the diamond-tipped head, the finer the grain. Keep this in mind according to the treatment to be performed.
- In addition to the grain, there are heads of different sizes. Select the appropriate measurement according to the area to be treated.
 - Body, especially for knees, elbows and stretch marks.


- Face and large spots

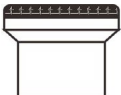- Scars, fine lines and imperfections

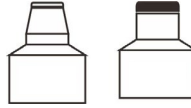- Eyes, nose, wrinkles, blackheads and acne.

CLEANING AND MAINTENANCE

- Before cleaning the device, unplug it from the mains.
- Clean the device with a soft, dry cloth. Do not use liquids to clean it.
- After each use, the various components and accessories must be sterilised and left to dry.
- For maintenance or repairs of any part of the product, please contact your distributor or specialised technician.

RECYCLING

If you need to dispose of the product or any of its parts or accessories, please check and respect the laws and regulations applicable to waste management and environmental protection.

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