vibbel

0009351

Khona Press

Pressotherapy / Presoterapia / Pressothérapie / Pressoterapia / Pressotherapie

EN	2
ES	6
FR	10
IT	14
DE	18

DDUUEETT EUROPE would like to thank you for choosing Vibbel products. We hope that this instruction manual will help you get the most out of this product.

If you have any doubts regarding its use and/or operation, please contact your distributor.

BEFORE USE

- Please read the entire instruction manual carefully before using the product. Take note of any warnings included in this manual.
- Always keep this instruction manual handy and if you give the device to another person, hand over the manual as it is an essential part of the product.
- Some components of the product may differ from the image that appears in catalogues, websites, brochures, etc.
- In order to enhance the product, the manufacturer reserves the right to modify its features and/or components without prior notice.

CORRECT USE OF THE DEVICE

- This product has been designed to be used by beauty professionals at beauty salons and should not be used for other purposes.
- Please note that any misuse of the device entails the loss of the guarantee.
- Do not make any modification to the device without the authorisation of DDUUEETT EUROPE SL.
- DDUUEETT EUROPE SL shall not be held liable for any material or personal damage due to the misuse of the product.

CONTRAINDICATIONS

- This product should not be used on:
 - People who suffer from deep vein thrombosis.
 - People with pacemakers or acute heart failure.
 - People with high blood pressure.
 - Pregnant women.
 - People with wounds or injuries in the area to be treated.
 - People with pain, irritation, pins and needles, cramp or numbness in the area to be treated.
- Not recommended for cancer patients or people with an overactive thyroid.
- The following should check with their doctor:
 - People with skin disorders.
 - People undergoing medical treatment or a diagnostic process.

PRECAUTIONS

- Before using it for the first time and after each use, check the device and the mains cable for any damage.
- Unplug the device from the mains before starting to clean it or undertaking any type of maintenance.
- Do not use the product outdoors or in places where it may be affected by humidity or inclement weather.
- If the power supply cable is defective, unplug the appliance immediately and inform your dealer or specialised technician.
- Position the power supply cable properly so that people cannot trip over it and make sure you plug it into an earthed socket.
- At the end of the working day, disconnect the device and unplug the power supply cable from the mains.
- If you detect any anomaly in the operation of the product, please contact your distributor or specialised technician.

SAFETY INDICATIONS

- The device and accessories should be kept out of the reach of children.
- Only plug the device into a properly installed socket. The local mains voltage must correspond to the technical specifications of the device.
- The power socket must be accessible in case you have to quickly unplug the device at any given time.
- Always pull out the actual plug and not the mains cable.
- Never tamper with any electrical component or open up/repair any component of the device.
- Only your distributor or specialised technician should repair the device.
- Do not switch on the device if there is visible damage to it or to the mains cable.
- Never touch the device or the mains cable with wet hands.

CONTENTS

- Main body x 1 unit
- Power cable x 1 unit
- Suit: 2 boots, 2 legs, 2 arms and waist
- Mat and air pressure tubes
- Fuse x 2 units
- Instruction manual x 1 unit

TECHNICAL SPECIFICATIONS


Product	Pressotherapy
Reference	0009351
Voltage	AC 200-230 V~
Input frequency	50 Hz
Power	50 W
Dimensions	44 x 31 x 15 cm

COMPONENTS

- Main body

- Suit

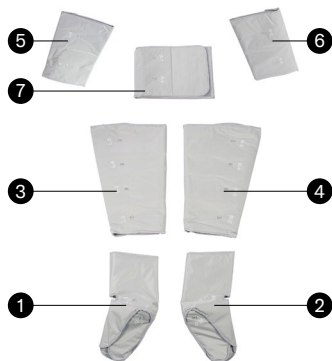- 1 Right boot
- 2 Left boot
- 3 Right leg
- 4 Left leg
- 5 Right arm
- 6 Left arm
- 7 Waist


- Mat and air pressure tubes

CONTROL PANEL


- 1 On / Off button
- 2 Display screen and parameter control

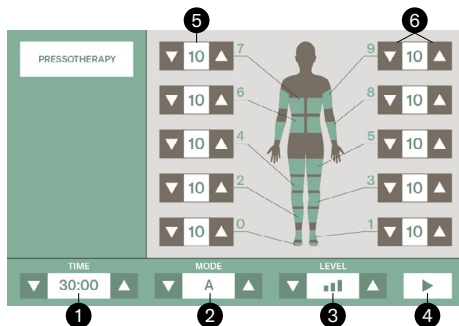
- 1 Time display and adjustment
- 2 Mode display and adjustment
- 3 Pressure intensity display and adjustment
- 4 Start or Pause the operation of the machine
- 5 By pressing the square until “-” appears, the corresponding area will not be inflated
- 6 Inflation time adjustment for each area (in seconds)

- Mode C: Foot to arm sequence, with all sections kept inflated and followed by the simultaneous inflation of all the sections.

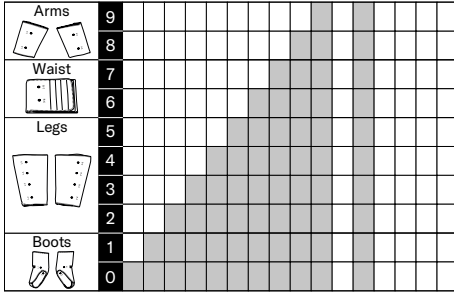- Mode D: Sequence with simultaneous activation of legs and arms.

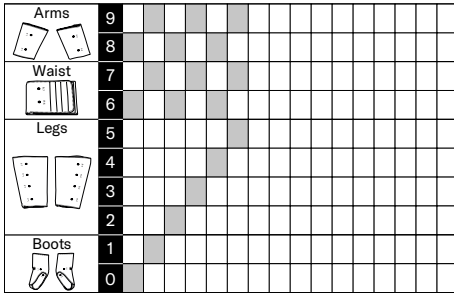CLEANING AND MAINTENANCE

- Before cleaning the device, unplug it from the mains.
- Clean the device with a soft, dry cloth. Do not use liquids to clean it.
- After each use, the various components and accessories must be sterilised and left to dry.
- For maintenance or repairs of any part of the product, please contact your distributor or specialised technician.

RECYCLING

If you need to dispose of the product or any of its parts or accessories, please check and respect the laws and regulations applicable to waste management and environmental protection.

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